

# karma

A Better Wheelchair, A Better Fit

**KARMA MEDICAL PRODUCTS CO., LTD.**

NO.2363,Sec.2,University Rd.,Min-Hsiung Shiang,Chia-Yi 621,Taiwan

**T +886-5-206 6688 F +886-5-206 7788**

Email: [globalsales@karma.com.tw](mailto:globalsales@karma.com.tw)

[www.KarmaMedical.com](http://www.KarmaMedical.com)



# karma

A Better Wheelchair, A Better Fit



Power Wheelchair

## eflexx

*Travel Freely and Lightly*

The lightweight power wheelchair is easy to assemble and light to travel with. Moreover, its high quality frame not only ensures the user's safety but also offers the customer the ultimate maneuverability.



CE  ISO 13485

Karma reserves the right to modify information herein without further notice 111501000106 Rev. 04 Dec. 2021





eFlexx has been created to provide the user with the best possible experience available.

Simple, Light and Efficient

Selectable battery range: large capacity/small capacity

- Standard 26.1Ah (657.7Wh) lithium battery: provides an extended driving range.
- Optional 11.5Ah (290Wh) lithium battery: travel friendly, adheres to IATA\* rules.

No need to take out the battery, eFlexx can be easily folded into a compact form for ease of transportation.

\*Each passenger can carry one spare battery with a capacity of less than 300 Wh or two less than 160 Wh. Please note that airlines might have their own rules and regulations with regard to flight safety, please follow their instructions accordingly.

Battery charge has never been this easy



Battery installation requires just one simple step.

Lightweight & removable battery allows you to charge it anywhere and anytime.

Benefit of Lithium Battery

Ecological friendly, highly efficient and long lasting. The battery weight (3.9 kg) is only 1/6 the weight of the traditional lead-acid batteries.



Safety First

Reinforced double cross-braces

An engineering design that gives you a comfortable seating experience and ensures the user's safety.



Transit and Felf-propel

Two versions are available : we have different rear wheel sizes can accommodate different make it the same line



EFL-20  
20" aluminum rear wheels

Switchable manual/power mode



EFL-14  
14" rear wheels  
Easy to carry



Flexible and Adjustable Design

Configurable seat width and seat depth, swing-away footrests, backrest angle, armrest height are all adjustable. You can fully personalize your wheelchair according to your needs.



Angle Adjustable Backrest  
-4°/0°/4°/8°/12°/16°

Easy to operate angle adjustment from -4° to 16°, giving users a more comfortable experience and flexible choices.



Universal Design

The bright color of the key lockers are self-explanatory; users can easily operate the functions without any hesitation.



Swing-away Controller System

Swing-away controller system makes it easier to approach a desk-top. The new safety feature prevents hands and fingers from pinching accidents.



Foldable Push Handles

The new foldable push handles allow for minimal effort in folding and pushing the wheelchair, as well as giving eFlexx a sleek and streamline look.



Swing-away Footrests

Improves convenience in tight areas and allows easier access to a chair or a bed.



Flip-back Armrests with Integrated locking system

Allows for easier transfer in and out of the chair.



Height Adjustable Armrests

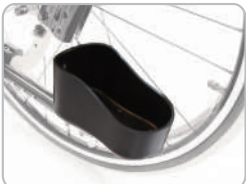
Tool-free adjustable armrest height (21-29cm)

OPTIONS & ACCESSORIES

Karma flexx's optional accessories are compatible with eFlexx.



USB Charger



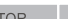



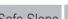




Cane Holder



Axle-traveling Elevating Legrest

SPECIFICATIONS

Product Name	Seat Width	Seat Depth	Seat Height	Backrest Height	Seat Angle (X°)	Backrest Angle (X°)	Armrest Height	Legrest Length	Overall Length	Overall Width	Overall Height	Folded Length*Width*Height
EFL-14	36/41/46	41.5/44/46.5	48.5	41	7°	6-26°	21-29	32.5-42.5	101.5-106.5	56/61.5/66.5	86	81.5/86.5/91.5*50.5*73
EFL-20	36/41/46	41.5/44/46.5	48.5	41	7°	6-26°	21-29	32.5-42.5	112.5-117.5	59/64.5/69	86	83.5/89.5/94.5*52*71.5

Product Name	 MOTOR W	 BATTERY		Front / Rear Wheels	 Max Safe Slope	 KM/HR	 RANGE KM	 Max. Obstacle Height	 Turning Radius	 MAX	Weight w/o Accessories	Weight with Standard Batteries	Weight w/o Batteries
EFL-14	250W	26.1Ah (657.7Wh) 11.5Ah (290Wh)	Linx	8" / 14"	6°	8	20(26.1 Ah) 10(11.5 Ah)	3.5	94.5	100	27	35	30.9
EFL-20	250W	26.1Ah (657.7Wh) 11.5Ah (290Wh)	Linx	8" / 20"	6°	8	20(26.1 Ah) 10(11.5 Ah)	3.5	94.5	100	29	37.8	33.7

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

Available Color Silver Black