

karma



- ✓ Lightweight
- ✓ Ergonomic Back and Seat
- ✓ Lockable Attendant Brakes


Ergo 3

The Next Generation of
Lightweight Mobility

Karma's new Ergo Series wheelchairs feature the innovative 'Double S-Ergo' seat and backrest system, which offers unparalleled ergonomic support. This patented design ensures a stable, upright posture, keeping users comfortably aligned with their backs straight and heads held high.



Easy Line of Sight

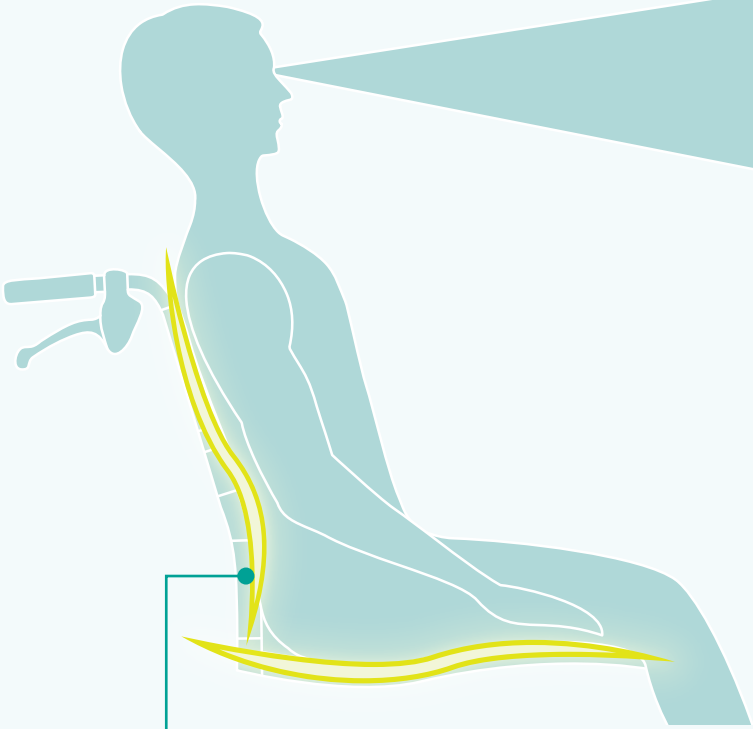

Senior wheelchair users often have reduced upper body control, which can lead to a slouched posture or a downward-facing head position.

Balanced posture starts with the right angles—**The S-Ergo Back System** combines scientifically validated seating angles with targeted lumbar support to:

- Promote a stable, healthy posture
- Enhance the user’s field of vision
- Improve social engagement



Ergo Lumbar Support



The S-Ergo Back Inserts are designed to provide targeted lumbar support on lightweight wheelchairs. They can be adjusted and expanded based on the user height and weight.




World-Renowned S-Ergo Seat

Follows the body’s anatomic curves to evenly distribute weight from the hip bone to the lower thigh, maximizing comfort and encouraging proper posture.

 **TRADITIONAL**

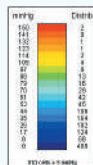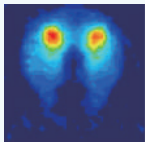

 **S-Ergo**

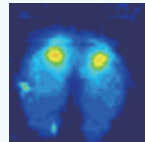

Intuitive Design

Designed with both users and caregivers in mind, all movable parts are highlighted in high-visibility orange for easy identification and operation.



Push & Pull Brake System



Flip-Up Armrests



Detachable Swing-Away Footrests



Multi-Functional Side Guards

Discreet, easily accessible pockets are elegantly integrated into the inner sides of the side guards for added convenience, while the outer sides feature reflective accents to enhance visibility and safety outdoors.



Drum Brake System (Optional)

The drum brake design delivers stable and powerful braking force, avoids direct tire wear to extend tire life, and provides an easy operating experience.



Easy Fit Legrest System (Optional)

Quick, tool-free height adjustment for comfort and support. Easy use and caregiver-friendly.

Lightweight Innovation

Despite its many innovative features, Ergo 3 is one of the lightest wheelchairs on the market, with a mass of the heaviest part just 9 kg (about 20 pounds) in its base configuration.



Comfort in Style



Color Options



Accessories

- Ergo Anti-Tippers
- Basic Anti-Tippers
- Pelvic Belt
- Cane Holder
- Oxygen Tank Holder
- IV Pole
- Foldable Headrest



Ergo 3 Specifications

Seat Width (mm)	410					460					510				
Seat Depth (mm)	435					435					435				
Front Caster Diameter (in)	6		7			6		7			6		7		
Rear Wheel Diameter (in)	14	20	14	20	24	14	20	14	20	24	14	20	14	20	24
Seat Surface Height at Front Edge (mm)	435		485			435		485			435		485		
Backrest Height (mm)	435					435					435				
Footrest to Seat Distance (mm)	385 - 460					385 - 460					385 - 460				
Seat Plane Angle (°)	6					6					6				
Seat to Backrest Angle (°)	104					104					104				
Max. User Weight. (kg)	125*					125*					125*				
Min. Transport Weight (kg)	Q 8.0 F 10.1	Q 8.2 F 11.3	Q 8.4 F 10.5	Q 8.4 F 11.5	Q 8.4 F 12	Q 8.3 F 10.4	Q 8.5 F 11.6	Q 8.7 F 10.8	Q 8.7 F 11.8	Q 8.7 F 12.3	Q 8.6 F 10.7	Q 8.8 F 11.9	Q 9 F 11.1	Q 9 F 12.1	Q 9 F 12.6
Min. Total Mass (kg)	11.9	13.1	12.3	13.3	13.8	12.2	13.4	12.6	13.6	14.1	12.5	13.7	12.9	13.9	14.4
Min. Overall L * W * H (mm)	1050 x 575/595 x 860 1050 x 575/595 x 910					1050 x 625/645 x 860 1050 x 625/645 x 910					1050 x 675/695 x 860 1050 x 675/695 x 910				
Min. Folded L * W * H (mm)	Q 770 x 245 x 660/F14 770x 285 x 700 F20 780x 310 x 700/F24 820x 310 x 700					Q 770 x 245 x 700/F14 770x 285 x 735 F20 780x 310 x 735/F24 820x 310 x 735					Q 770 x 245 x 740/F14 770x 285 x 780 F20 780x 310 x 780/F24 820x 310 x 780				

*The Ergo 3 is rated for users up to 125 kg. However, when equipped with Push & Pull attendant brakes, the maximum user weight is reduced to 115 kg. Karma reserves the right to modify information herein without further notice. "/ means optional : "-" means adjustable ± 1 cm ± 0.5kg Rev.01 Dec. 2025



Ergo 3

KARMA MEDICAL PRODUCTS CO., LTD.
T +886-5-206 6688
Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com