

karma



- ✓ Ergo Back
- ✓ Lightweight
- ✓ More Weight Capacity

Ergo 3

The Next Generation of Lightweight Ergonomic Mobility

At its core, Karma Ergo Series combines ergonomic excellence and lightweight mobility. Now it takes a leap forward by creating a natural balance that spans from head to shoulders, down to the lower back. This clever design frees users from the effort of maintaining an upright head position, enhancing their social life with more joy and vitality.


At the heart of this transformation lie the patented Ergo Back and market-proven S-Ergo Seat. Discover the future of seating with the Ergo 3 - where form, function, and well-being converge in perfect harmony.



Advanced Ergo Seating System

Elderly wheelchair users tend to have weaker upper body control which leads to slouched posture or downward-facing head position. The advanced Ergo Seating System addresses these common issues by incorporating scientifically proven seating angles.

- Effectively promoting balance in body posture
- Facilitating an expanded field of vision
- Elevating social interactions
- Mitigating chronic sitting discomfort

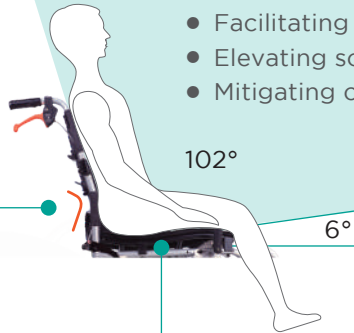

Patent Ergo Back

A set of curved and adjustable inserts provide direct lumbar support and maintains the spine's natural alignment, for maximum comfort and support in a lightweight wheelchair.



Market Proven S-Ergo Seat

Matches the anatomic curves of the body, to evenly distribute users' weight from the hip bone to the lower thigh to maximize comfort and encourage proper posture.

Push & Pull Brake

Just push to part and pull to unlock; No need to bend down to park the wheelchair!



Side Guard Pocket

A dedicated space to keep your personal items close and secure, with reflective material to maximize safety on the road.



Easy Transfer Armrest and Legrest

The flip-up armrests and detachable swing away footrest enable easy transfers on and off the wheelchair.



Scan to learn more about Ergo 3

- Accessories
- Ergo Anti-tipper
 - Additional Ergo Back Insert
 - Pelvic Belt
 - Anti-Tippers
 - Cane Holder
 - IV Pole Holder
 - Oxygen Tank Holder

Frame Color



Pearl Silver Silver Black Rose Red

Upholstery Color



Gray Dark Gray Purple Orange Yellow

Ergo 3 Specifications

Seat Width* Seat Depth (inch)	16*17 / 18*17 / 20*17
F/R wheels (inch)	7 / 14 、 7 / 24
Seat Surface Height at Front Edge (mm)	485
Backrest Height (mm)	435
Armrest Height (mm)	220
Calf Length (mm)	335 - 435
Seat to Backrest Angle (X°)	102
Max. User Weight (kg)	115 (14F/Q) 125 (24F/Q)
Total Mass (kg)	11.4 - 14.2
Overall L x W x H (mm)	1080 x 570 x 915 1080 x 700 x 915
Folded L x W x H (mm)	800 x 310 x 740 800 x 330 x 760

Remark : Karma reserves the right to modify information herein without further notice. " / " means optional : " - " means adjustable $\pm 1\text{ cm}$ $\pm 0.5\text{ kg}$
Rev.00 Aug.2023



KARMA MEDICAL PRODUCTS CO., LTD.

T +886-5-206 6688

F +886-5-206 7788

Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com