

karma



- ✓ Lightweight
- ✓ Ergo Seat
- ✓ Lockable Attendant Brakes

Ergo Lite 2 TL



The Same Lightweight Freedom, Now with a Better Fit

The Ergo Lite 2 TL takes our all-time best-selling model to the next level, offering expanded size options, a higher weight capacity, and improved seat-to-floor height. It's the perfect blend of comfort, practicality, and portability for a wider range of users.



Everyday Users



Caregivers




Travelers

Who is **Ergo Lite** series for?



Ergo Series Design Concept

Matches the anatomic curves of the body, to evenly distribute users' weight from the hip bone to the lower thigh to maximize comfort and encourage proper posture.

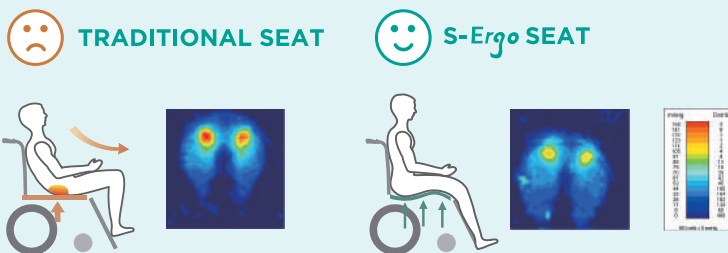

Detachable, Swingaway Footrests

Allows for easy transfer, on and off the wheelchair and detaches with no tools required.



Push & Pull Brake

Karma's Push & Pull brake system allows attendants to safely park the wheelchair without bending down to reach the wheel locks.



Easy to Transport and Store

Convenient and suitable for elderly caregivers to take care of elderly users. Easy to carry out and fold for storage.



Crash Tested

The Ergo Lite 2 TL has been crash tested in accordance with ISO 7176-19, making it suitable for use as a seat in a motor vehicle.

With reflective material to maximize safety on the road.

Accessories

- Foldable Headrest
- Pelvic Belt
- Harness
- Anti-Tippers
- Cup Holder



Tool-free quick-release rear wheels

Specifications

±1 cm ± 0.5kg

Seat Width*Depth (inch)	16 x 17	18 x 17	20 x 17
F/R wheels (inch)	7 / 14		7 / 22
Seat Surface Height at Front Edge (mm)	510		
Backrest Height (mm)	435		
Calf Length (mm)	410 - 460		
Seat to Backrest Angle (X°)	98		
Max. User Weight (kg)	115		
Total Mass (kg)	9.8 - 10.4		11.1 - 11.9
Overall L x W x H (mm)	min. 990-1005 x 565 x 950 max. 1015-1030 x 665 x 950	min. 990-1005 x 575 x 950 max. 1015-1030 x 700 x 950	
Folded L x W x H (mm)	min. 647 x 240 x 745 max. 720 x 270 x 785	min. 720 x 260 x 690 max. 755 x 280 x 785	
Crash Test	✓		

The actual product specifications may differ from the data in the table by ±1 cm ± 0.5 kg.
Remark : Karma reserves the right to modify information herein without further notice.
"/" means optional : "-" means adjustable Rev.00 Apr. 2025

