

karma A Better Wheelchair, A Better Fit



karma

A Better Wheelchair, A Better Fit

KARMA MEDICAL PRODUCTS CO., LTD.

NO.2363 Sec,2 University Rd.,Min-Hsiung Shiang,Chia-Yi 621,Taiwan

T +886-5-206 6688 F +886-5-206 7788

Email: globalsales@karma.com.tw
www.KarmaMedical.com




Power Wheelchair



CE   ISO 9001

Ergo Nimble
KP-10.3S

Independence with Ergo comfort

Ergo nimble is an affordable rear wheel drive power chair. It features the patented S-Ergo Seating System and promotes a comfortable seating and a good posture.

Compact design makes it ideal for confined spaces mobility and short distance journey.



Mainly for indoor use and short distance journeys



The 75cm turning radius for easy entry and exit of small elevators, making indoor activity effortless.



Tool-free disassembly. The maximum weight per piece is 26kg, can be easily placed in a car trunk, making outdoor travel possible and simple.



Flip-back armrest



Adjustable armrest height



Tool-free adjustable backrest angle (95°-109°)

Optional Accessories



The optional quick release steering controller enables maneuvering switches from the user to the attendants in seconds.



Swing-away and detachable footrest



Adjustable armrest height

Equipped with captain seat for extra riding comfort




KP-10.2 Captain Seat Version



Tool-free disassembly, making outdoor trip possible and simple


TRADITIONAL SEAT

Causes high risk of decubitus ulcers because of pressure points concentrating on hipbone area.

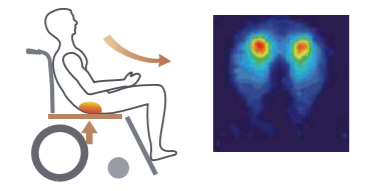


S-Ergo SEAT

The S-Ergo system (intelligent s-shaped ergonomic seating) provides efficient pressure relief by spreading weight over a greater area. This also increases stabilization and reduces sliding.

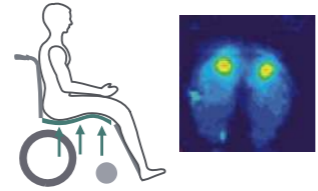

HOW DOES IT WORK?

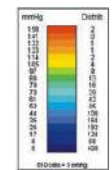

Designed according to the body curve, S-Ergo Seating is capable of relieving pressure, stabilizing seating position and reducing downward sliding.

Specification

44	54	103	62	94	38-46	manual: 5°, 19°	Normal: 200W Max: 420W	22AhX2	200X50/12"	VR2	75	8°	5	6km/hr	28km	115kg	61.5kg	37.5kg

Karma reserves the right to modify information herein without further notice.

Unit : cm ±1cm ±0.5kg