


Ergo Recline

Reclining Transfer Wheelchair

The Ergo Recline offers a near-horizontal recline with flip-back detachable armrests for convenient lateral transfers and excellent comfort.

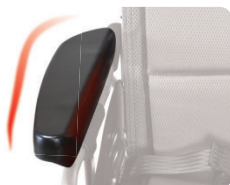

Boat Armpad

Wider than traditional armpads and featuring a comfortable, contoured shape, they are easily re-positioned to accommodate self-propulsion and transport.



Ergonomic Headrest with 3D Pillow

Provides lateral and occipital support to relieve neck and shoulder pressure without restricting vision. Velcro backing allows you to adjust the headrest position to create the perfect fit.

Patented

Multi-functional Armrests

Flip-up, height-adjustable, and fully detachable armrests accommodate a wide range of use cases, such as changing clothes and transferring.



Elevating Legrests

Angle-adjustable, with calf pads, for improved comfort and support.



Hydraulic Pressure Support System

The Ergo Recline's recline system contains two pneumatic gas struts, which provide a smooth, controlled tilting motion and actively assist the attendant return the backrest to the upright position.

With a 161° reclining angle, the seating pressure can be redistributed across the entire body.





SPECIFICATION

Nominal seat width (mm)*	415 , 465
Nominal seat depth (mm)**	420
Front caster / Rear wheel size (in)	8 x 14 / 8 x 20
Front seat height (mm)	500
Backrest height (mm)	495
Armrest to seat distance (mm)	220 - 300
Footrest to seat distance (mm)	365 - 440
Max. user weight (kg)	115
Mass of the heaviest part (kg)	13.3 / 13.4 / 14.2 / 14.3
Total mass (kg)	21.1 / 21.2 / 21.8 / 21.9
Overall L x W x H (mm)	Min. 1190 x 590 x 700 / Max. 1815 x 690 x 1300
Overall folded L x W x H (mm)	790 x 340 x 940 / 850 x 400 x 940



Scan to learn more
about Ergo Recline

* Nominal seat width is measured the width of seat cushion at 120mm in front of the backrest.
This value is mainly used in the product ordering process.

** Nominal seat depth is measured from the front end of the seat cushion to the back cushion.
This value is mainly used in the product ordering process.

- The actual product specifications may differ from the data in the table by $\pm 1\text{ cm}$ $\pm 0.5\text{ kg}$.
- The actual product specifications may vary according to different configurations.
- Karma reserves the right to modify information herein without further notice.

Rev.04 Jun . 2025

KARMA MEDICAL PRODUCTS CO., LTD.
T +886-5-206 6688
Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com