

Ergo Series

Two ergonomic models with our Double S-Ergo seat and back system—designed for maximum comfort and tailored to different user needs.



PATENTED

S-Ergo Seating System

Contours to the body's natural curves to evenly distribute pressure across the seat surface and reduce pressure points.




PATENTED

S-Ergo Back System

Incorporates clinically validated seat and back angles with targeted lumbar support to improve posture, stability, and overall comfort.

Easy Line of Sight

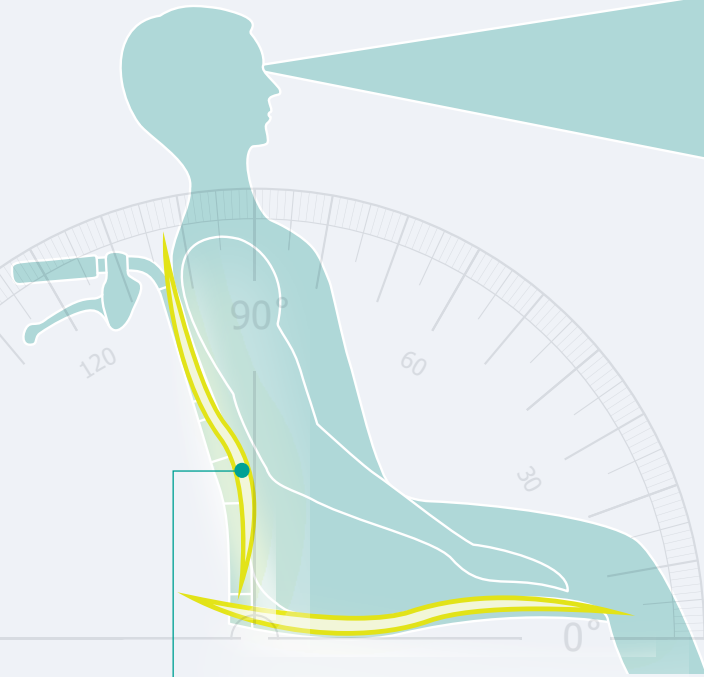

Senior wheelchair users often have reduced upper body control, which can lead to a slouched posture or a downward-facing head position.

Balanced posture starts with the right angles—**The S-Ergo Back System** combines clinically validated seating angles with targeted lumbar support to:

- Promote a stable, healthy posture
- Enhance the user's field of vision
- Improve social engagement

Ergo Lumbar Support

The S-Ergo Back Inserts provide targeted lumbar support with a customizable fit that adjusts to the user's height and weight, enhancing comfort and posture. This advanced feature is available exclusively on the Ergo 2 and Ergo 3.



Intuitive Design

Designed with both users and caregivers in mind, all movable parts are highlighted in high-visibility orange for easy identification and operation.



Attendant Brakes



Flip-Up Armrests



Swingaway Detachable Footrests



Ergo Anti-Tippers

The Next Generation of Lightweight Mobility



Karma's new Ergo Series wheelchairs feature the innovative 'Double S-Ergo' seat and backrest system, keeping users comfortably aligned with their backs straight and heads held high.

Ergo 3

Flip-Up armrests and detachable, swing-away footrests for easy side transfers.



Ergo 3 Self-propel



Ergo 3 Transit

Ergo 2

Fixed armrests and detachable, swingaway footrests make this an ideal option for forward transfers.



Multi-Functional Side Guards

Discreet, easily accessible pockets are integrated into the inner sides of the side guards, while the outer sides feature reflective accents to enhance safety outdoors.



Lightweight Innovation

Despite its many innovative features, Ergo 3 is one of the lightest wheelchairs on the market, with its heaviest part weighing only 8 kg in its base configuration.

Accessories

- Ergo Anti-Tippers
- Basic Anti-Tippers
- Pelvic Belt
- Cane Holder
- Oxygen Tank Holder
- IV Pole
- Foldable Headrest



Tool-free Adjustable Footplates (Optional)



Drum Brake System (Optional)



Scan to learn more about Ergo 3

Series	Product Name	Nominal seat width (mm)	Nominal seat depth (mm)	Front caster / Rear wheel size (in)	Seat surface height at front edge (mm)	Backrest height (mm)	Armrest to seat distance (mm)	Footrest to seat distance (mm)	Backrest angle (X°)	Seat plane angle (X°)	
Ergo	Ergo 3	410	435	6 / 14 7 / 14	435 485	435	220	385 - 460	20	6	
		460									
		510									
		410									
		460									
		510									
		410									
		460									
	510										
	Ergo 2	Ergo 2	410	435	6 / 14 7 / 14	435 485	435	220	385 - 460	20	6
			460								
			510								
			410								
			460								
			510								
			410								
460											
510											

Leg to seat angle (X°)	Maximum user weight (kg)	Mass of the heaviest part (kg)	Total mass (kg)	Overall length with legrest (mm)	Overall width (mm)	Overall height (mm)	Folded length (mm)	Folded width (mm)	Folded height (mm)
105	115	8.2 / 8.6 (14Q) 10.3 / 10.7 (14F)	12.1 / 12.5	1010 - 1030	575	885 / 935	725 (Q) 725 (F)	245 (Q) 285 (F)	640/660 (Q) 665/715 (F)
		8.5 / 8.9 (14Q) 10.6 / 11.0 (14F)	12.4 / 12.8		625				690/710 (Q) 665/715 (F)
		8.8 / 9.2 (14Q) 10.9 / 11.3 (14F)	12.7 / 13.1		675				740/760 (Q) 720/760 (F)
	115 / 125	8.4 / 8.6 (20Q) 11.5 / 11.7 (20F)	13.3 / 13.5		595		725 (Q) 780 (F)	245 (Q) 310 (F)	640/660 (Q) 665/715 (F)
		8.7 / 8.9 (20Q) 11.8 / 12.0 (20F)	13.6 / 13.8		645				690/710 (Q) 665/715 (F)
		9.0 / 9.2 (20Q) 12.1 / 12.3 (20F)	13.9 / 14.1		695		740/760 (Q) 720/760 (F)		
		8.6 (24Q) 12.2 (24F)	14.0		595	725 (Q) 825 (F)	660 (Q) 715 (F)		
	8.9 (24Q) 12.5 (24F)	14.3	645		935	725 (Q) 825 (F)	710 (Q) 715 (F)		
	9.2 (24Q) 12.8 (24F)	14.6	695				760 (Q) 760 (F)		
	105	115	8.0 / 8.4 (14Q) 10.1 / 10.5 (14F)		11.8 / 12.2	1010 - 1030	575	885 / 935	725 (Q) 725 (F)
8.3 / 8.7 (14Q) 10.4 / 10.8 (14F)			12.1 / 12.5	625	690/710 (Q) 665/715 (F)				
8.6 / 9.0 (14Q) 10.7 / 11.1 (14F)			12.4 / 12.8	675	740/760 (Q) 720/760 (F)				
115 / 125		8.2 / 8.4 (20Q) 11.3 / 11.5 (20F)	13.0 / 13.2	595	725 (Q) 780 (F)		245 (Q) 310 (F)		640/660 (Q) 665/715 (F)
		8.5 / 8.7 (20Q) 11.6 / 11.8 (20F)	13.3 / 13.5	645					690/710 (Q) 665/715 (F)
		8.8 / 9.0 (20Q) 11.9 / 12.1 (20F)	13.6 / 13.8	695	740/760 (Q) 720/760 (F)				
		8.4 (24Q) 12.0 (24F)	13.7	595	935			725 (Q) 825 (F)	660 (Q) 715 (F)
8.7 (24Q) 12.3 (24F)		14.0	645	710 (Q) 715 (F)					
9.0 (24Q) 12.6 (24F)		14.3	695	760 (Q) 760 (F)					

Rev.04 Jun . 2025

KARMA MEDICAL PRODUCTS CO., LTD.
T +886-5-206 6688
Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com