



KARMA MEDICAL PRODUCTS CO., LTD.

NO.2363 Sec,2 University Rd.,Min-Hsiung Shiang,Chia-Yi 621,Taiwan

T +886-5-206 6688 F +886-5-206 7788

Email: globalsales@karma.com.tw
www.KarmaMedical.com



Karma reserves the right to modify information herein without further notice 111501000040 Rev. 00 Jan 2016

karma A Better Wheelchair, A Better Fit



Power Wheelchair

EVO LECTUS  

A Chair That Evolves with You

EVO LECTUS represents an evolution from years of experience in wheelchair assessments.

It's born to suit the exact specification of each market, and can be customized to fit the delicate needs of every single client.



EVO LECTUS[®]



LECTUS seat Full Scale Adjustable

The full-scale adjustable system that can be set up to the individual needs, including fully multi adjustable seat, backrest, armrests and legrest.



EVO chassis Compact and powerful

Equipped with a 4-wheel superb full suspension system which guarantees powerful and smooth driving performances in all terrains. With 45cm low seat height, 62cm in narrow total width, 88cm in length after flipping up the footplates, above all provide an agile using experience.



MAIN FEATURES

40 cm Power Elevating

Allow seat height from 45 cm up to 85cm.

180 Degree Power Recline + 180 Degree Power Legrest

Total flat position allow the user able to lay down straight.




50 Degree Power Tilting

Dynamic CoG design allow tilt up to 50 degree.

Modular Design

- Every function can be modular upgrade
- Multi-adjustable Armrest (5 axis) & Footrest (3 axis)

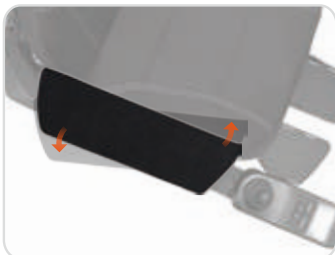

Modular Power Setting

Including dynamic COG tilt, biomechanical recline and legrest, seat lift.

Advanced Electronics

EVO LECTUS is compatible with R-net to provide various possibilities.



OPTIONAL ACCESSORIES



Parallel Swing Away-down



3D Adjustable Headrest



Hip Support



Lateral Support




Calf Support



Transfer Bar

- Push Handle
- Retractable Pelvic Belt
- LED Light

Available Colors

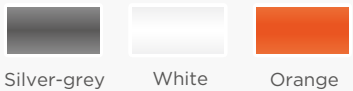

SPECIFICATION

Seat Width	Seat Height	Seat Depth	Overall Length	Overall Width	Overall Height	TILT	WHEELS (F/R)	MOTOR	BATTERY	WHEELS (F/R)	R-Net
42-52cm	45cm	40-60cm	120cm	62cm	112cm	0° - 50°	4 Wheel	320W	75Ah	14" / 9"	R-Net
Turning Radius	Max Safe Slope	KM / HR	RANGE	MAX	Weight with Batteries	Weight w/o Batteries	±1cm ±0.5kg -Karma reserves the right to modify information herein without further notice. -For more details, please refer to the Specification Pages No. 35-38.				
75cm	10°	7cm	6Km 10Km	40Km	140kg	165kg					