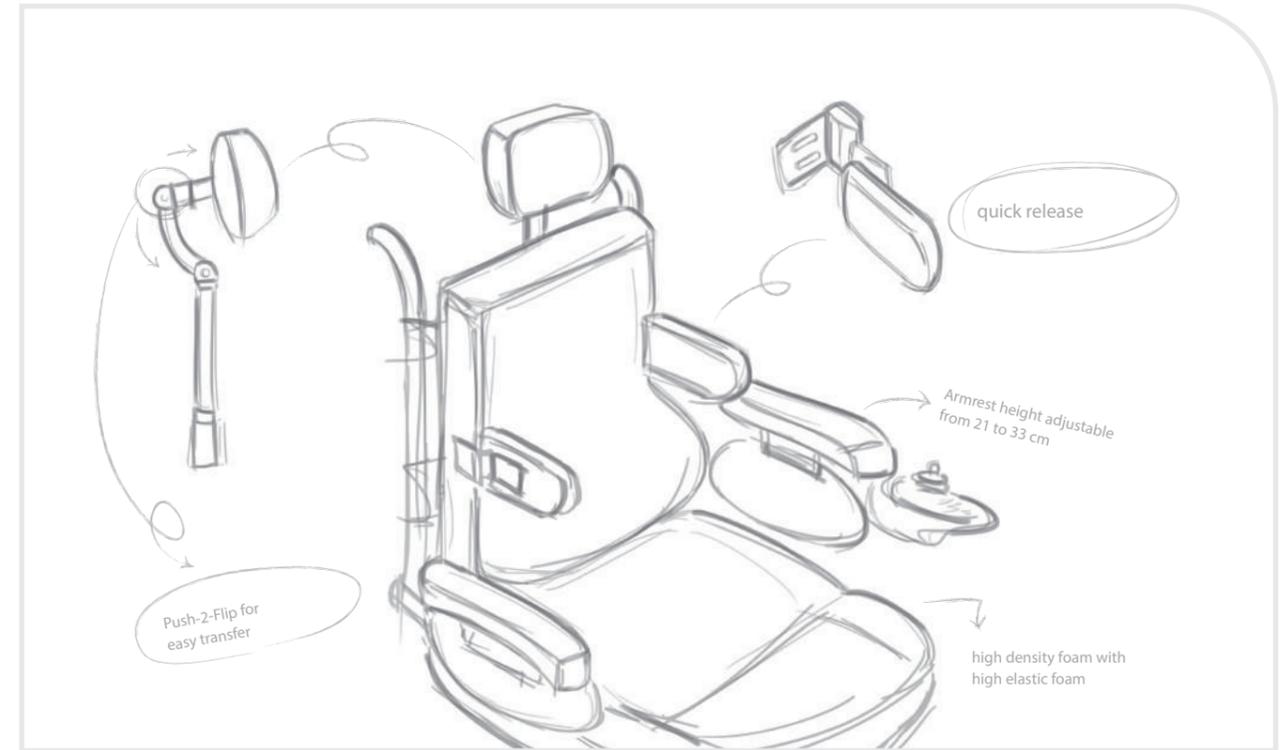


karma

 A Better Wheelchair, A Better Fit

karma
A Better Wheelchair, A Better Fit

KARMA MEDICAL PRODUCTS CO., LTD.

NO.2363 Sec,2 University Rd.,Min-Hsiung Shiang,Chia-Yi 621,Taiwan

T +886-5-206 6688 F +886-5-206 7788

Email: globalsales@karma.com.tw
www.KarmaMedical.com




Seating System

Kiss

 Karma Integrated Seating System

Born To Fit The Most Challenging Demands / Born To Fit Better Solutions

At Karma, every part of the users' lifestyle is taken into account while the seating system is being created. The brand new **KISS** (Karma Integrated Seating System) offers numerous thoughtful functions that make daily life easier. The unique Push-2-Flip design is built in all supporting accessories to lighten the load of care givers while transferring.

A fully adjustable seat frame, multiple power features and various accessories have made KISS the best choice in assisted mobility devices.



A New Chapter of the Balance of Performance and Individual Clinical Needs

- Multiple power features- 50° tilt, 170° shearing force reduction recline, 30cm lift, power elevating legrest.
- Modular supporting accessories for advanced positioning.
- Push-2-Flip design for easy operation and transfer.
- Easily adapted to different chassis.



Power Lift

-Conquer the challenges of different height in daily life by 30cm lift.



Power Tilt

- 50° tilt provides perfect share of pressure loading from bottom to back and enhances the posture stability.
-The advanced center of gravity system features the most safety and stability.



Power Recline

-With power elevating legrest and 170° shearing force reduction recline, KISS provides the best seating comfort and pressure relief.




Fully adjustable seating system that grows with the user

KISS allows for width, depth, back height and angle adjustments to be made easily and quickly in the field, which can grow with the client.

	KISS	KISS Junior
Seat Width	16"- 20"	12"-16"
Seat Depth	16"- 20"	12"-16"
Backrest Height	20"- 22"	16"-18"



Backrest Pan

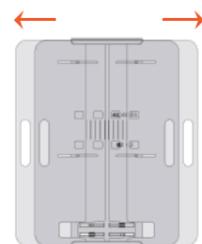


Seat Pan




Pull-2-slide Seat Pan



Sliding Rail

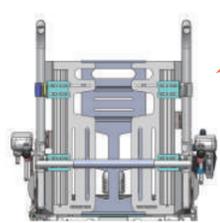

Backrest Height

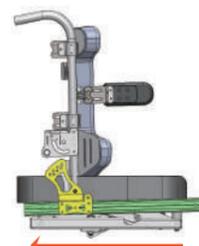

Seat Depth

Adaptable for Leon and MGN power base



Features & Options

Headrest



Flip back and detachable headrest



Curve adjustable headrest

- Easy to adjust by just one tool.
- Height and multiple-direction adjustable to reduce the pressure of head and shoulder.
- Push-2-Flip for easy transfer.

Swing-away Controller



Perfect combination of user friendly mechanism allows closely access to the environment.

Upholstery Fabric




Waterproof, breathable and stretchable upholstery for pressure relief.

Options of Armrest



Height adjustable flip-back armrest

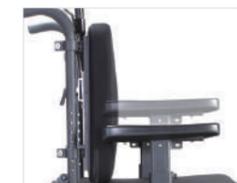

T-shape height adjustable & detachable armrest



Reclining armrest



Ergo armpad offers comfortable support for elbow.



Contour armpad provides optimal support of weak upper limbs.



Padded armpad offers comfort and protection to elbows and wrists.


Options of Legrest



Vertical legrest



Standard legrest

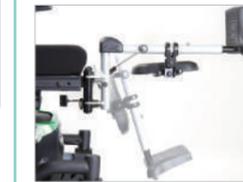

Elevating legrest



Power elevating legrest



Angle adjustable footplate



Multi-adjustable Footplate

Modular Supporting Accessories

- Modular supporting accessories with more than a hundred combinations for advanced positioning.
- Sliding rail for easy adaption and adjustment.



Trunk support



Hip support



Pommel



Push-2-Flip design

With one simple button, care givers can easily swing away supporting accessories and complete the transfer.