


KM-BT10

Strong, Compact, and Adjustable

KM-BT10 is a fully configurable bariatric wheelchair; the foldable dual cross brace and seat pan allow for a weight capacity of 250 kg and a maximum seat width of 760 mm, all while retaining the portability and practicality of a standard wheelchair. KM-BT10 is also highly adjustable to create and maintain the perfect fit.

Prescription based seat and backrest system.

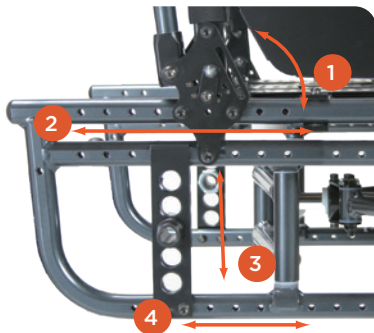

- 1 Back angle
- 2 Seat depth
- 3 Seat to floor height
(Available with different wheel and caster package)
- 4 Center of gravity



The foldable aluminum dual cross brace frame make the KM-BT10 one of the lightest and most compact bariatric wheelchairs on the market.



Solid Seat Pan


Provides rigidity and support, to prevent the Hammock Effect and ensure a smooth ride.



Rigidifying Push Bar

Strengthens the wheelchair frame and provides attendants with an ergonomic surface to comfortably and efficiently maneuver the wheelchair.

Multi-directional Adjustable Footplates

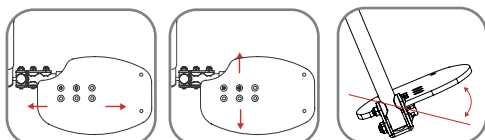

SPECIFICATION

Nominal seat width (mm)*	560 / 610 / 660 / 710 / 760
Nominal seat depth (mm)**	460 / 510 / 560 / 610
Front caster / Rear wheel size (in)	7/8 x 24Q
Front seat height (mm)	390 -500
Backrest height (mm)	410 - 510
Armrest to seat distance (mm)	250
Footrest to seat distance (mm)	320 - 370
Max. user weight (kg)	250
Mass of the heaviest part (kg)	22.9
Total mass (kg)	31
Overall L x W x H (mm)	Min. 1180 x 780 x 800 / Max. 1230 x 980 x 1010
Overall folded L x W x H (mm)	Min. 900 x 440 x 800 / Max. 950 x 440 x 1010

* Nominal seat width is measured the width of seat cushion at 120mm in front of the backrest.
This value is mainly used in the product ordering process.

** Nominal seat depth is measured from the front end of the seat cushion to the back cushion.
This value is mainly used in the product ordering process.

- The actual product specifications may differ from the data in the table by ± 1 cm ± 0.5 kg.
- The actual product specifications may vary according to different configurations.
- Karma reserves the right to modify information herein without further notice.

Rev.04 Jun . 2025