



Manual Wheelchair

SOMA 215

Compact, Easy Transit, and Lightweight Solution

SOMA 215 is designed for easy day-to-day transportation. The patented swivel legrests can be swung away with one hand in one second, providing a safer frontward transfer. The chair is also very compact that it can be fitted into a car trunk as small as 72 x 29 x 68cm, which is less than two carry-on suitcases.



SOMA 215


Patented Swivel Legrest

Lift and Swing, just a simple action to free up leg space for transfer. The legrest is not detachable, so you will never lose it.



Detachable Seating System

Innovated Detachable Seating System makes it extremely easy to detach the upholstery so that you can clean it as often as you want.

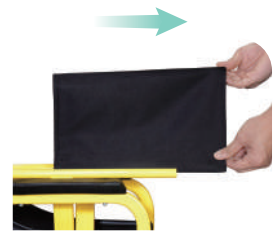

Compact and Lightweight

Its 9.9kg ultra-lightweight makes your day to day transportation effortless. And, the compact size is easy to fit into a small car trunk.



Seamless Elliptical Frame

The seamless elliptical-frame design reduces the amount of the welding points, improving the durability and the look of the chair.



Ergonomic Push Handle

The handle is designed according to the shape of hands. It is comfortable to hold and easy to push. The attendant brake can be added on upon request.



Anti-Slippery Footplate

An anti-slippery footplate provides you a stable and comfortable base for your foot.



Energy-saving Tipping Lever

Energy-saving tipping lever, allows caregivers to overcome obstacles easier.



14" and 20" Rear Wheel Options













SPECIFICATION

Upholstery Color



Frame Color



Model	F/R Wheels											Max. User Weight
SOMA 215	6"x14"	40 / 45	39	60 / 60.5	26	22	93	45	39	86	9.9kg / 10 kg	100kg
	6"x20"	40 / 45	39	60 / 60.5	29	22	95	45	39	86	10.9kg / 10.9kg	100kg

±1cm ±0.5kg -Karma reserves the right to modify information herein without further notice. 111501000130 Rev.01 Mar. 2018



KARMA MEDICAL PRODUCTS CO., LTD.

T +886-5-206 6688

F +886-5-206 7788

Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com