


SOMA 215

Compact and Lightweight

At only 10.2 kg, the SOMA 215 is perfect for daily use. Never lose a footrest again; the patented mechanism swings completely out of the way without detaching, making transfers, storage, and transportation easy and worry-free.

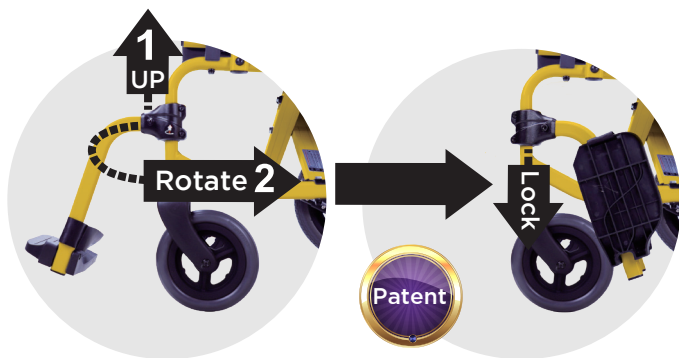

Never Lose a Footplate Again!

The patented swingaway footrests are easy to swing away without the need of detaching it, making everyday transfers, storage, and transportation smooth and worry-free.



Seamless Elliptical Frame

Crafted using minimal components and welds, resulting in both a refined aesthetic and exceptional durability.



Detachable Upholstery

Intuitive design makes it easy to detach and clean the upholstery.



Ergonomic Push Handles

Offers a stable and comfortable operating experience for attendants.



Attendant Tipping Lever

Allow caregivers to overcome obstacles with ease.



Various sizes and colors to suit
your needs and style!



Optimized Seat Heights

The SOMA 215 features an accessible 43 cm seat-to-floor height, ideal for users under 162 cm or those who prefer to foot propel, while the SOMA 215 Tall's 48 cm seat height is better suited for users over 162 cm.



Self-propel

Model	Seat Width x Depth (mm)	Seat to Floor Height (mm)
SOMA 215	400 x 395	430
Self-propel 20"	425 x 395	
Transport 14"		
SOMA 215 TL	460 x 420	480
Self-propel 22"		
Self-propel 20"		
Transport 14"		



Transport



Anti-Slip Footplates

Improved grip allows users to better maintain proper foot position.

FRAME COLOR

- Bright Silver White
- Yellow
- Flat Black



Scan to learn more
about SOMA 215 TL



Series	Product Name	Nominal seat width (mm/in)	Nominal seat depth (mm/in)	Front caster / Rear wheel size (in)	Seat surface height at front edge (mm)	Backrest height (mm)	Armrest to seat distance (mm)	Footrest to seat distance (mm)	Backrest angle (X°)	Seat plane angle (X°)
SOMA	SOMA 215	425	405	6 x 14	430	400	220	345 - 400	11	8
		425		6 x 20						6
		460	430	7 x 14	480	430		430 - 474	10	5
		460		7 x 20						
		460		7 x 22Q						

Leg to seat angle (X°)	Maximum user weight (kg)	Mass of the heaviest part (kg)	Total mass (kg)	Overall length with legrest (mm)	Overall width (mm)	Overall height (mm)	Folded length (mm)	Folded width (mm)	Folded height (mm)
95	100	10.2	10.2	940	605	855	630	300	690
		11.2	11.2	950			715	310	
110		10.5	10.9	1010	640	940	680	255	760
		11.7	12	1020			775	280	
		8.5	12.4	1045	660		785	300	

Rev.04 Jun . 2025

KARMA MEDICAL PRODUCTS CO., LTD.
T +886-5-206 6688
Email : globalsales@karma.com.tw

karma
www.KarmaMedical.com