


**karma** A Better Wheelchair, A Better Fit

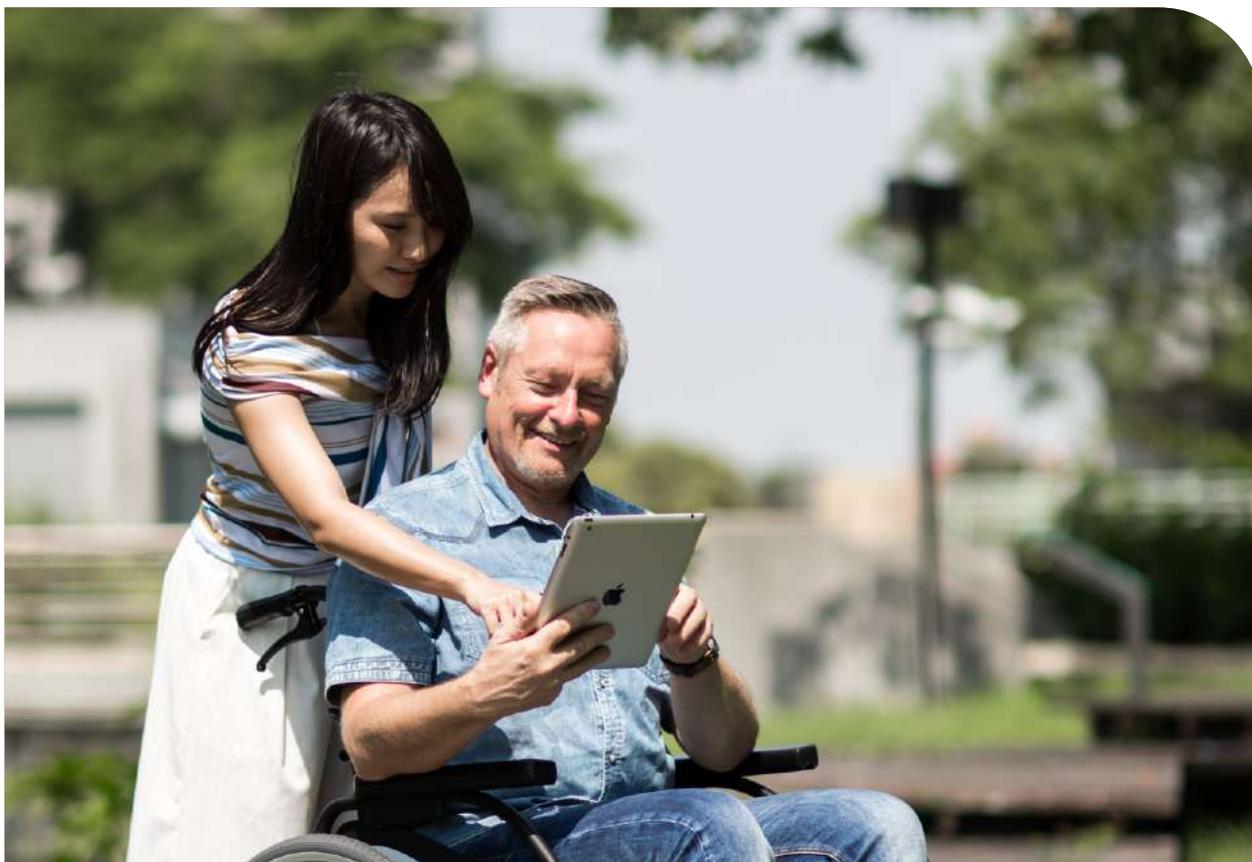

Manual Wheelchair

**★ STAR 2**

***Stylish, Lightweight,  
and Easy to Transfer***


Karma STAR 2 is designed to change the image of traditional wheelchair. Besides its catching design, it provides excellent maneuverability and numerous ergonomic details. It also makes front transfer safer by offering patented swivel footrest.





## Increase front transfer safety by creating larger leg space in one second

Did you know it is very likely for the user to fall when moving to and from a wheelchair? The moment that the user is transitioning, is when the most attention to safety is needed. The Karma STAR 2 offers the patented swivel legrests, which free up leg space, for both attendant and the user, and increase your everyday safety.

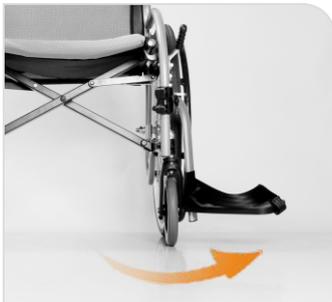


### ◀ Patented Swivel legrest

Lift and Swing, just a simple action to free up leg space for transfer. The legrest is not detachable, so you will never lose it.



### ◀ Detachable Seating System

Innovated Detachable Seating System makes it extremely easy to detach the upholstery, so that you can clean it as often as you want. It also enhanced the support on the surface, which can avoid the "hammock effect."

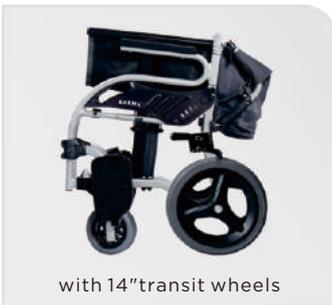

with 14" transit wheels

### ◀ Compact, Easy to Carry

Foldable backrest and the swivel-away footrest reduce the overall size of the chair which makes it easier to store and to carry around.



- The seamless elliptical-frame design reduces the amount of the welding points, improving the durability and the look of the chair.
- Energy-saving tipping lever, allows caregivers to overcome obstacles easier.
- Anti-skid protected footplate, preventing the chair from damaging furniture.

## SPECIFICATION

width	F/R wheels											Max. User Weight
16"	7"x24"	40	43	62	29	22	100	50	42	92	13.1 kg	115 kg
18"	7"x24"	45	43	67	29	22	100	50	42	92	13.5 kg	115 kg
16"	7"x 14"	40	43	60	27	22	94	50	42	92	12.5 kg	115 kg
18"	7"x 14"	45	43	65	27	22	94	50	42	92	12.9 kg	115 kg

111501000125 Rev.00 Dec.2017

±1cm ±0.5kg -Karma reserves the right to modify information herein without further notice.



KARMA MEDICAL PRODUCTS CO., LTD.

T +886-5-206 6688

F +886-5-206 7788

Email : globalsales@karma.com.tw

**karma**  
www.KarmaMedical.com