

# VIP 2



## Curved Headset

The standard headrest is fully adjustable, including the height, horizontal position, depth, angle, and curvature, ensuring a perfect fit.




## J-Hooks

VIP 2 features ergonomic rotating J-shaped push handles, which allow attendants to operate and maneuver the wheelchair without reaching too far, even when the wheelchair is fully tilted and reclined.

## Detachable Footrests

The standard detachable, swingaway footrests can be quickly removed or repositioned to clear the space in front of the wheelchair, allowing users to propel with their feet or perform forward transfers.

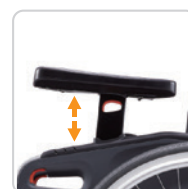

## Ergonomic Armrest

Height-adjustable and detachable for comfort and convenience.

## OPTIONS & ACCESSORIES

### Biomechanical Elevating Legrests

The VIP 2's elevating legrests automatically extend during elevation, allowing the user's legs to stretch naturally, reducing pressure behind the knees.

- Tension Adjustable Seat
- Angle Adjustable Footplates
- Quick-Release Axles
- Pneumatic Tires
- Adjustable Curved Headrest
- Chest Belt (Velcro Type)
- Chest Belt (Plastic Buckle)
- Pelvic Belt
- Hemi Armpad
- Wooden Tray Table

- Amputee Supports
- Vertical Footrests
- Backrest Bag
- Net Bag
- Cane Holder



Scan to learn more about VIP 2

### Weight Shifting With Comfort

Karma VIP 2 is designed to avoid the feeling of anxiousness that happens when going in and out of tilt. The unique weight shifting system allows the center of gravity to travel between front and rear wheels while tilting, providing the steadiest balance possible and a great sense of security for the user.

### Anti-Sliding Reclining System

VIP 2 minimizes sliding and shear forces experienced by the user by rotating along the user's hip joint, keeping them comfortable and secure through the full range of recline and tilt motions.



### Tension Adjustable Backrest

The internal structure of the backrest is formed by a series of individually-adjustable Velcro straps, which can create a user-specific fit to maximize comfort and support.



### Foldable and Durable

VIP 2 is light, foldable, and easy to fit into the trunk of a car. The rigidifying push bar strengthens the wheelchair frame and provides attendants with an ergonomic surface to more comfortably and efficiently maneuver the wheelchair.



### Modularized Design, Easy to Upgrade

VIP 2 is also compatible with most third-party accessories and backrests.



Series	Product Name	Nominal seat width (mm)	Nominal seat depth (mm)	Front caster / Rear wheel size (in)	Seat surface height at front edge (mm)	Backrest height (mm)	Armrest to seat distance (mm)	Footrest to seat distance (mm)	Backrest angle (X°)	Seat plane angle (X°)						
Positioning	VIP 2	400	415	6/7 x 14	470 / 490	510 / 560	220 - 280	320 - 420	13 - 43	1 - 36						
		400	465													
		450	415													
		450	465													
		400	415	6 / 7 x 20							470 / 490	510 / 560	220 - 280	320 - 420	13 - 43	1 - 36
		400	465													
		450	415													
		450	465													

Leg to seat angle (X°)	Maximum user weight (kg)	Mass of the heaviest part (kg)	Total mass (kg)	Overall length with legrest (mm)	Overall width (mm)	Overall height (mm)	Folded length (mm)	Folded width (mm)	Folded height (mm)
100	115	17.8	24.4	1070	630	1380	775	355	680
		18.0	24.7	1120			825		
		18.1	24.9	1070	680		775		
		18.3	25.2	1120			825		
		15.4(20Q) 19.0(20F)	25.4	1070	660		840	380	715
		15.6(20Q) 19.2(20F)	25.6	1120			890		
		15.6(20Q) 19.2(20F)	25.6	1070	710		840		
		15.8(20Q) 19.4(20F)	25.8	1120			890		

Rev.04 Jun . 2025