

# VIP 515



## Tilt-in-Space with Foldable Frame

The VIP 515 is one of the few foldable tilt-in-space wheelchairs on the market, making it an ideal choice for users who have difficulty maintaining posture but do not want to be limited by their wheelchair.



## Adjustable Flip-up Armrests

A proper armrest height is helpful to relieve loads on the neck, shoulders, arms, and back. No tools required. The armrest also flips up for easier transfers on and off the wheelchair.




## Pneumatic Pressure Support System

The VIP 515's tilt-in-space system contains two hydraulic gas struts, which provide a smooth, controlled tilting motion and actively assist the attendant return the wheelchair to 0°.

Tilt range: 0° - 35°

## Foldable and Compact

The VIP 515 has an aircraft-grade aluminum frame, which is lightweight and foldable without using tools, offering greater convenience for storage and travel.

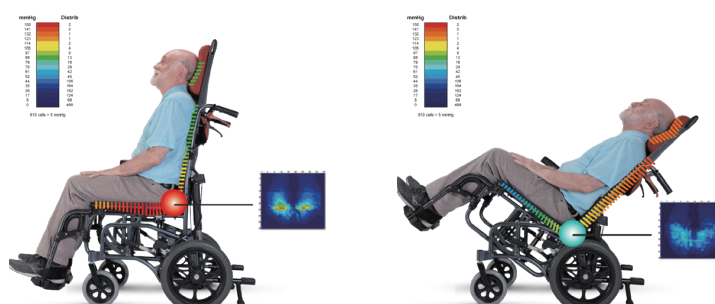

## Tilt-in-Space

Easily reposition the user by tilting the chair to redistribute pressure over a larger surface area, reducing the risk of developing pressure sores and back pain.



## Easy Transfer

Flip-up armrests and swing away footrests improve convenience and provide easier access to the chair.

## SPECIFICATION

Nominal seat width (mm)*	405 , 455
Nominal seat depth (mm)**	415
Front caster / Rear wheel size (in)	7 x 14 / 7 x 20
Front seat height (mm)	490 , 485
Backrest height (mm)	460
Armrest to seat distance (mm)	240 - 320
Footrest to seat distance (mm)	405 - 480
Max. user weight (kg)	115
Mass of the heaviest part (kg)	15.1 / 15.3 / 16.2 / 16.4
Total mass (kg)	18.9 / 19.3 / 20.0 / 20.3
Overall L x W x H (mm)	Min. 1070 x 635 x 955 / Max. 1070 x 710 x 1330
Overall folded L x W x H (mm)	Min. 725 x 350 x 760 / Max. 790 x 375 x 760

\* Nominal seat width is measured the width of seat cushion at 120mm in front of the backrest. This value is mainly used in the product ordering process.

\*\* Nominal seat depth is measured from the front end of the seat cushion to the back cushion. This value is mainly used in the product ordering process.

- The actual product specifications may differ from the data in the table by  $\pm 1\text{ cm} \pm 0.5\text{ kg}$ .

- The actual product specifications may vary according to different configurations.

- Karma reserves the right to modify information herein without further notice.



Scan to learn more  
about VIP 515

Rev.04 Jun . 2025

KARMA MEDICAL PRODUCTS CO., LTD.  
T +886-5-206 6688  
Email : globalsales@karma.com.tw

**karma**  
www.KarmaMedical.com